

INTEGO. PERFORMANCE YOU CAN RELY ON.



SIRONA.COM

The Dental Company

sirona.

YOUR RELIABLE PARTNER

INTEGO is the best treatment center of its class. It is based on the innovative power and award-winning design of our existing treatment centers. As a reliable partner for an efficient practice, INTEGO has all the features that provide optimal support: ergonomic design for better workflows, high quality standards in processing and design, and a great value – this is why INTEGO is our performance class. **Enjoy every day. With Sirona.**

ERGONOMIC DESIGN

Designed for intuitive handling and optimized, integrated workflows. More on pages 06 – 09

QUALITY MADE IN GERMANY

Defined quality assurance processes, durable materials, and sturdy, future-proof design. More on pages 18 | 19

GREAT VALUE

Excellent performance and innovative functions at a great price, customized options for great value. More on pages 04 | 05

CONTENTS

- 04 Versions and operating concepts
- 06 Ergonomic design
- 10 Easy operation
- 12 Flexible, future-proof design
- 14 Patient comfort
- 16 Hygiene functions
- 18 Quality by Sirona
- 20 Individual configuration

GREAT VALUE YOU CAN RELY ON.

Your work day is complex and places many different demands on the treatment center you choose. INTEGRO gives you the perfect solution: Choose from two models that you can configure flexibly. Whether you select an INTEGRO or INTEGRO pro, a whip arm or a hanging hoses design – every option offers optimal performance at a fair price.



INTEGRO

- Treatment center with a compact water unit for optimal use of space.
- EasyPad user interface for easy, intuitive operation.
- Compact assistant element with good freedom of movement.
- Controlled water quality with a self-contained fresh water supply bottle.
- Available as the TS model with hanging hoses or CS with whip arms.

INTEGRO pro

- Integrated, comprehensive hygiene concept for optimal results.
- EasyTouch touchscreen and integrated foot control concept for easy, intuitive operation.
- Flexible assistant element suitable for two- or four-handed treatment.
- Numerous equipment options with a great value.
- Available as the TS model with hanging hoses or CS with whip arms.

pro

This icon on the following pages of the brochure indicates equipment that is available only with the **INTEGRO pro**.

ERGONOMIC DESIGN YOU CAN RELY ON.

The treatment center is the heart of your practice. INTEG0 perfectly combines appearance and functionality – so you can work comfortably under optimal conditions. In addition to their harmonious design, coordinated components and variety of color combinations, INTEG0 treatment centers meet the strict criteria of our ergonomics concept. INTEG0 fits the style of your practice – and your taste.

THE COLOR COMBINATIONS

THE ELEGANT WORLD

Classic, timeless combinations for an especially high-class look.

Carbon

Aubergine

Platinum

Pacific

THE NATURAL WORLD

Soft, warm shades for a soothing treatment atmosphere.

Basalt

Rosewood

Cherry

Marble

THE VITAL WORLD

Vivid colors and bright tones for a modern, lively atmosphere.

Sapphire

Spring

Orchid

Lagoon

AWARD-WINNING DESIGN

Sirona has earned prestigious design prizes for many of its products, including the red dot design award, the iF Product Design Award Gold, the iF communication design award, and the GOOD DESIGN Award.

Explore these color choices and design for yourself with the color configurator at [SIRONA.COM/INTEGO](https://www.sirona.com/intego).



PERFECT RESULTS.

Every treatment is a complex procedure. Our ergonomic products provide optimal access to the patient and smooth treatment workflows with a high level of patient comfort. INTEGRO combines every dimension of ergonomics to make treatment easier for you. Achieve perfect results even at the end of a long working day – with your expertise and the support of our treatment center.



1	INTUITIVE SITTING	2	COMFORTABLE POSITIONING	3	OPTIMUM VISIBILITY	4	INTEGRATED WORKFLOW	PERFECT RESULTS
	More freedom of movement Thin backrest for ideal freedom of leg movement, flexible positioning of the foot control	Easy adjustment of the treatment center Flexible height adjustment, swiveling cuspidor	Perfect illumination of the treatment area LEDview or LEDlight operating light, LED instrument light	Intuitive operation EasyPad user interface or EasyTouch touchscreen	<ul style="list-style-type: none">■ Clinically perfect■ Esthetically perfect■ Perfect for the success of your practice		Integrated tabletop device functions Integrated torque control for endodontics, ApexLocator on the EasyTouch	
	Optimized instrument placement Perfect coordination of dentist and assistant elements, flexible positioning of the additional tray in the whip arm and hanging hoses model.	Best spine support Ergonomic, comfortable upholstery. Synchronized backrest and seat movements for optimal patient positioning. No headrest readjustment needed when the backrest is tilted.	Flexible headrests Flat or double articulating headrest	Integrated tabletop device functions Integrated torque control for endodontics, ApexLocator on the EasyTouch			Modern equipment Innovative instruments, modern patient communication system SIVISION, USB, SD card and Ethernet interface	
	Comfortable working stools HUGO, CARL, PAUL	Greater patient comfort Gentle chair lift drive, armrests, Vario footrest	Perfect intraoral images SiroCam F or SiroCam AF, 22" HD monitor at the treatment center					

INTUITIVE OPERATION

Difficult treatments require your full concentration. This is why INTEGRO is intuitive – with two new, easy-to-understand user interfaces and innovative instruments for stress-free treatment. So you can concentrate on what's really important – your patient and a perfect treatment result.



EasyPad

Clear and intuitive user interface standard in both INTEGRO and INTEGRO pro.

BL ISO E micro-motor

Especially light, compact and with an ISO interface. Speed and torque are ideal for endodontic treatment.




Integrated torque control

The motor stops automatically when the set maximum torque is reached. This ensures greater safety and a perfect workflow for endodontic treatment.



EasyTouch

This user interface can be integrated optionally in the INTEGRO pro. More operating comfort thanks to intuitive touchscreen technology and more operating and display options.

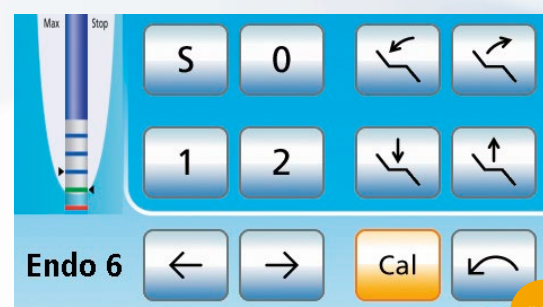

ApexLocator

An ApexLocator can also be integrated into the EasyTouch to allow permanent control of the working depth during endodontic treatment.

SPRAYVIT E
Optimal visibility with LED light.

SIROSONIC L SCALER
Handpiece and tips can be sterilized; controlled, straight oscillation. With electronic power and water control.

SIROCAM F OR AF
Easy-to-install plug & play cameras for detailed intraoral images. The SIROCAM AF with autofocus is especially easy to operate.

FLEXIBLE INTO THE FUTURE


In a good practice, customized settings add to your comfort. This is why you can configure the INTEG0 exactly the way you need for optimal support of your treatment workflow. But flexibility means even more to us – and you can prepare INTEG0 with technical options and interfaces so you are ready for future applications.




MANY MODELS

CHOICE OF TWO HEADRESTS

The flat headrest can be adjusted quickly and flexibly. Optionally, you can select a double articulating headrest that is easier to position and allows you optimal visibility even in difficult-to-see areas.

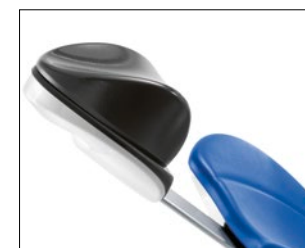

Flat headrest

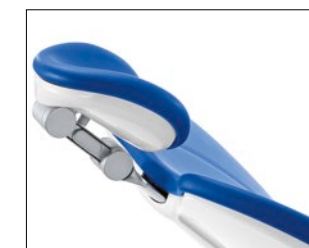

Double articulating headrest

TWO OPERATING LIGHT OPTIONS – ALWAYS OPTIMAL ILLUMINATION

Two choices of LED lights for optimal illumination of the operating field. They can be positioned flexibly and operated hygienically using a no-touch sensor.



LEDlight



LEDview

ASSISTANT ELEMENT: COMFORT OR COMPACT

Both assistant elements are very space-saving and ensure maximum freedom of movement. With the Comfort model, you can work comfortably even without an assistant.



Comfort assistant element



Compact assistant element

Even more variety.
See page 23 for more configuration options to adapt the INTEG0 optimally to your needs.

FUTURE-PROOF DESIGN



Integrated USB interface

INTERFACES: FLEXIBLE INTO THE FUTURE

A network interface in the patient chair and a USB port in the dentist or assistant element allow you to integrate new functions and network your practice – troubleshooting by your service technician is also much easier.



SIVISION 22" AC monitor

SIVISION 22" AC monitor

With this monitor, you can explain intraoral and X-ray images, software and planning views, videos and PowerPoint presentations to your patient directly at the treatment center.

IMPRESSIVE PATIENT COMFORT

Your patients will also appreciate this comfort: The INTEG0 patient chair has an ergonomic design and, together with the headrest, ensures that your patients are especially comfortable. The chair lowers far enough for even shorter patients to can get in and out easily. And there is no doubt – patients who are comfortable feel more secure and remain loyal to your practice.



DOUBLE ARTICULATING HEADREST
Quickly adjustable to various treatment situations and patients with different statures.

PREMIUM UPHOLSTERY
The high-quality upholstery ensures comfort even during long treatments.

SYNCHRONIZED BACKREST AND SEAT MOVEMENTS FOR OPTIMAL PATIENT POSITIONING
The tilting of the backrest matches the movement of the patient's spine to prevent stretching and compression effects.

ARMRESTS
Your patients can rest their arms comfortably. The left armrest is used for support when getting in and out; the right one is moved out of the way.

LIFTING SYSTEM
The unique chair lift drive gently moves patients into the proper position.

VARIO FOOTREST
The Vario footrest that comes with the Premium upholstery can be easily extended to provide even very tall patients with ergonomic positioning.



SAFE HYGIENE

Dentists all over the world have the highest demands for the safety of their patients and practice team. The INTEGRO hygiene concept is designed so you can comply quickly and easily with all hygiene regulations. The assistant element and water unit designs make hygiene procedures even easier and less time consuming. All of this guarantees perfect infection control – for you, your patients, and your team.

INTEGO

The compact water unit and assistant element can be adapted to many different demands and requirements by freely selectable suction systems and hygiene options.



Fresh water bottle
For controlled water quality. Can also be filled with DENTOSEPT solution for permanent disinfection and manual sanitization.



Automatic suction hose cleaning
Cleaning adaptor for automatic purging of saliva ejector and spray aspirator.




Smooth surfaces
The largely gap-free design of the surfaces ensures maximum infection control in the shortest time possible.



INTEGO pro

Comfort water unit and assistant element are equipped with a comprehensive hygiene concept. You achieve the highest hygiene standards quickly and easily.

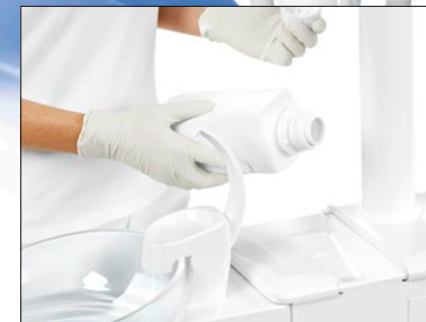

Integrated disinfection system
Permanent disinfection of all water lines with DENTOSEPT. Automatic sanitizing function for regular intensive disinfection. Can also be used with self-contained water supply.



Integrated cleaning adapters
For automatic purging of instrument hoses (AutoPurge) and sanitization of the water lines.



Chemical suction hose cleaning
Cleaning agent container for chemical purging of the suction lines. The concentration of cleaning agents can be adjusted individually.

QUALITY YOU CAN RELY ON.



Designed and produced in Germany. Using precisely defined processes that ensure greatest reliability.

We don't just want to satisfy our customers, we want to impress them. This is why quality, performance, and service are key for us. Sirona treatment units combine high quality with real innovations that perceptibly enhance treatment comfort and workflow.

INTEGO OFFERS HIGH ADDED VALUE AND TOP QUALITY – IN EVERY RESPECT

Reliability and durability

Precise, elaborate quality tests and robust materials make INTEG0 a reliable treatment center with an above-average service life.



High quality

INTEGO immediately shows our attention to detail. The materials used – robust, shapeable plastics, solid steel and aluminum components on the inside, high-quality glass and upholstery – meet the high quality demands of a premium product.



High standards

Sirona is a technology leader with more than 130 years of innovation history and sets high quality standards. Some 290 engineers and scientists in Research & Development not only pay close attention to quality; they also drive the innovation process.






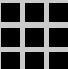




Future-proof design

At the Sirona innovation center, future trends are anticipated and products are prepared to be equipped for the technology demands of the future.



INDIVIDUAL CONFIGURATION

You can configure your INTEGO treatment center to fit your needs and meet your own desired great value.

BASIC UNIT	VERSION	USER PACKAGES	OPTION PACKAGE FOR DENTIST ELEMENT OR ASSISTANT ELEMENT
<div>INTEGO</div> <div>Patient chair<ul style="list-style-type: none">Premium upholsteryFlat headrestPneumatic instrument foot control</div> <div>Dentist element<ul style="list-style-type: none">5 instrument positionsEasyPad user interfaceSupply for 2 drill drives</div> <div>Assistant element<ul style="list-style-type: none">Compact version on a swivel arm directly mounted to the water unit4 handpiece positionsSaliva ejector and HVE</div> <div>Water unit<ul style="list-style-type: none">Compact designStandard wet suctionReady for the integration of a fresh water bottle systemCuspidor with manual swivel function</div>	<div>INTEGO TS With hanging hoses</div> <div></div> <div><ul style="list-style-type: none">The height of the dentist element can be adjusted and held in place with a pneumatic brakeErgonomic design with flexibly positionable dentist element and movable tray holderAmple instrument area on the dentist element and tray holder</div>	<div> START</div> <div><ul style="list-style-type: none">Standard 3-way syringeLight source for dentist element2 turbine hosesLEDlight operating light</div> <div> ADVANCE</div> <div><ul style="list-style-type: none">C+ electronic foot controlSPRAYVIT E in the dentist elementLight source in the dentist element1 turbine hose1 BL ISO E micro-motor with supply electronicsLEDlight operating light</div> <div> INNOVATION</div> <div><ul style="list-style-type: none">C+ electronic foot controlSPRAYVIT E in the dentist elementLight source for dentist1 turbine hose1 BL ISO E micro-motor with supply electronicsSIROSONIC L scalerLEDview operating light</div>	<div> SIVISION ready</div> <div><ul style="list-style-type: none">USB interfaceNetwork interface</div> <div> SIVISION Start Without monitor</div> <div><ul style="list-style-type: none">SIVISION readySIROCAM FMonitor mount on the light support column</div> <div> SIVISION Start With monitor</div> <div><ul style="list-style-type: none">SIVISION readySIROCAM FMonitor mount on the light support column22" AC monitor</div> <div> SIVISION Advance With monitor</div> <div><ul style="list-style-type: none">SIVISION readySIROCAM AFMonitor mount on the light support column22" AC monitor</div>

BASIC UNIT	VERSION	USER PACKAGES	OPTION PACKAGE FOR DENTIST ELEMENT OR ASSISTANT ELEMENT
<div>INTEGO pro</div> <div>Patient chair<ul style="list-style-type: none">Premium upholsteryPneumatic instrument foot controlFlat headrest4-way foot switch</div> <div>Dentist element<ul style="list-style-type: none">5 instrument positionsEasyPad user interfaceSupply for 2 drill drivesInstrument light sourceSupply electronics for up to 2 motors and scaler</div> <div>Assistant element<ul style="list-style-type: none">Comfort version, on a support arm with wide swivel range4 handpiece positionsSaliva ejector and HVELight source for SPRAYVIT E</div> <div>Water unit<ul style="list-style-type: none">Comfort designReady for the integration of an automatic disinfection deviceIntegrated cleaning adaptersCuspidor with manual swivel functionWet suction with multi-station valveIntegrated water flush suction hose cleaningWater heater</div>	<div>INTEGO pro TS With hanging hoses</div> <div></div> <div><ul style="list-style-type: none">The height of the dentist element can be adjusted and held in place with a pneumatic brakeErgonomic design with flexibly positionable dentist element and movable tray holderAmple instrument area on the dentist element and tray holder</div> <div>INTEGO pro CS With whip arms</div> <div></div> <div><ul style="list-style-type: none">Ideal for the widest range of treatment situations, whether you are working alone or in a teamOptimum range of instruments: up to 90 cmThe height of the dentist element can be adjusted and held in place with a pneumatic brake</div>	<div> START pro</div> <div><ul style="list-style-type: none">C+ electronic foot controlSPRAYVIT E in the dentist element1 turbine hose1 BL ISO E micro-motorLEDlight operating light</div> <div> ADVANCE pro</div> <div><ul style="list-style-type: none">C+ electronic foot controlSPRAYVIT E in the dentist element1 turbine hose1 BL ISO E micro-motorSIROSONIC L scalerEasyTouch user interfaceLEDview operating light</div> <div> INNOVATION pro</div> <div><ul style="list-style-type: none">C+ electronic foot controlSPRAYVIT E in the dentist elementThird drill drive1 turbine hose2 BL ISO E micro-motorsSIROSONIC L scalerEasyTouch user interfaceSPRAYVIT E in the assistant elementLEDview operating light</div> <div> ENDO pro</div> <div><ul style="list-style-type: none">Integrated torque control kitApexLocatorENDO 6:1 contra-angle handpiece</div>	<div> SIVISION ready</div> <div><ul style="list-style-type: none">USB interfaceNetwork interface</div> <div> SIVISION Start Without monitor</div> <div><ul style="list-style-type: none">SIVISION readySIROCAM FMonitor mount on the light support column</div> <div> SIVISION Start With monitor</div> <div><ul style="list-style-type: none">SIVISION readySIROCAM FMonitor mount on the light support column22" AC monitor</div> <div> SIVISION Advance With monitor</div> <div><ul style="list-style-type: none">SIVISION readySIROCAM AFMonitor mount on the light support column22" AC monitor</div>

YOUR INDIVIDUAL OPTIONS

Here you can view the complete list of all individual options, and a checklist of what is included in the basic units.

Options	INTEGO	INTEGO Pro
Patient chair		
2 treatment positions, S and O position	●	●
Number of user-specific settings	2	4
Central junction box	●	●
Chair lift system 140 kg	■	–
Chair lift system 185 kg	□	●
Flat headrest	■	■
Double articulating headrest	□	□
Premium upholstery	■	■
4-way foot switch	○	●
Booster cushion for small patients	○	○
Armrest, right	○	○
Armrest, left	○	○
Valve set for external devices	○	○
Pressure reduction kit	○	○
Cover installation kit	○	○
Separate junction box	○	○
Base plate for uneven floors, adapts for C-series installation to be replaced	○	○
Dentist element		
Height-adjustable support arm with pneumatic brake	●	●
5 instrument positions	●	●
Hydro-pneumatic supply for 2 drill drives w/o fiber optics	●	●
Pneumatic instrument foot control	■	■
C* electronic foot control	□A, I	□S, A, I
EasyPad user interface	●	■
EasyTouch user interface	–	□A, I
SPRAYVIT E with light	□A, I	□S, A, I
SPRAYVIT E with light (water on the left)	□A, I	□S, A, I
Standard 3-way syringe	□S	–
Turbine hose with light	○2xS, A, I	○S, A, I
Borden turbine hose	○	–
BL ISO E micro-motor with hose	○A, I	–
BL ISO E micro-motor with hose (ApexLocator ready)	–	○S, A, I
Second BL ISO E micro-motor (ApexLocator ready)	–	○I
SIROSONIC L scaler	○I	○A, I
Third drill drive	○	○I
Satelec Mini LED curing light	–	○
Intra-oral light source kit for turbine and SPRAYVIT E	○S, A, I	●
Torque control kit	○	○E
ApexLocator	–	○E
ENDO 6:1 contra-angle handpiece	○	○E
Intraoral X-ray viewer	○	○
Purge and manual sanitization function	○	–
Purge function	–	●
Additional TS tray holder for 2 standard size trays	○	○
Silicone pad for TS dentist element	○	○
Additional CS tray holder for 1 standard size tray	○	○
Additional CS tray holder for 2 standard size trays	○	○
Silicone pad for CS dentist element	○	○
Removable silicone sleeves for handles	○	○
Additional camera holder for TS dentist element	○	○

Options	INTEGO	INTEGO Pro
Assistant element		
Compact version, swivel mounted directly on the water unit	■	–
Comfort version, on a support arm with wide swivel range	□	●
4 handpiece positions	●	●
Saliva ejector and HVE	●	●
SPRAYVIT E with light	□	□I
SPRAYVIT E with light (water on the left)	□	□I
Standard 3-way syringe	□	–
2nd HVE kit	○	○
Satelec Mini LED curing light	○	○
Touch pad	○	○
Light source for SPRAYVIT E	○	●
Silicone pad for Comfort assistant element	○	○
Water unit		
Compact design, prepared for fresh water bottle	●	–
Comfort design, prepared for integrated disinfection device	–	●
Cuspidor with manual swivel function	●	●
Reduced price without cuspidor	□	–
Prepared for standard wet suction	■	–
Prepared for Standard wet suction w/automatic suction control	□	–
Multi-station valve	□	■
Dürr cuspidor valve	□	□
Automatic separator	□	□
Amalgam separator	□	□
Air jet pump (Air Venturi)	□	–
Fresh water bottle w/o switch to mains water supply	□	–
Fresh water bottle with switch to mains water supply	□	–
Integrated disinfection system with permanent waterline disinfection and automatic sanitization function	–	○
Water heater	○	●
Water flush suction hose cleaning system	○	■
Chemical suction hose cleaning system	–	□
Operating light options		
LEDlight	□S, A	□S
LEDview	□I	□A, I
Adaptation for HELIODENT Plus unit model	○	○
Patient communication		
USB interface in the dentist or assistant element	○SiR, SiS, SiA	○SiR, SiS, SiA
SIROCAM F	○SiS	○SiS
SIROCAM AF	○SiA	○SiA
Monitor mount on the light support column	○SiS, SiA	○SiS, SiA
22" AC monitor [with integrated power supply unit]	○SiS*, SiA	○SiS*, SiA
Network interface	○SiR, SiS, SiA	○SiR, SiS, SiA

- Included in the basic unit
- Option
- Variation, part of standard unit
- Optional variation
- Not possible/applicable

- SiR SIVISION ready package
- SiS SIVISION Start package with/without monitor
- SiS* SIVISION Start package with monitor
- SiA SIVISION Advance package
- S START/START pro package
- A ADVANCE/ADVANCE pro package
- I INNOVATION/INNOVATION pro package
- E ENDO pro package



HUGO working stool



CARL working stool

PAUL working stool
(pictured with additional foot ring)Additional swiveling tray holder for
CS dentist elementAdditional tray holder for TS dentist
element

Standard 3-way syringe



Satelec Mini LED curing light



X-ray viewer

Additional holder for camera
(for TS dentist element)Comfort assistant element with user
interface

C* electronic foot control



4-way foot switch



Valve set for external devices



Adaptation for HELIODENT unit model



22" AC monitor



INTEGO without cuspidor



LEDlight



LEDview



SIROCAM F



SIROCAM AF

ALWAYS AT THE FOREFRONT OF INNOVATION!

As global innovation leader for dental equipment, we continuously invest in research and thus in the future of modern dentistry. By networking digital technologies with integrated solutions and optimizing the treatment workflow, we can create improved treatment results, more comfort and safety for the patient as well as time and cost savings in everyday work. The combination of constant innovative power and globally growing sales and service structures makes Sirona the global market leader trusted by thousands of practices and labs around the world. **Enjoy every day. With Sirona.**



CAD/CAM systems

From pioneer to new standard. For almost 30 years we have been developing digital dentistry and creating new possibilities for the future practice and lab.



Imaging systems

Best image quality with the lowest dose. More than 100 years of developing x-rays for the dental practice make us the number 1 innovation partner.



Treatment centers

The business card of modern practices. We are striving to create the ideal ergonomic and innovative center. Individually tailored to the well-being and demands of the patient and dentist.



Instruments

Advantages that speak for themselves. We make sure that we provide the right balance of proven quality, individual ergonomics and innovative technology for user-friendly work.



Hygiene systems

Competence that gives you safety. When it comes to hygiene in the practice, we do not take any shortcuts.